



## Product Information Guide for COVID-19 Vaccines and Associated Products

**Contents**

Product Information Guide Summary .....2

Centrally Distributed Vaccines and Ancillary Kits.....3

    Vaccines (Centrally Distributed) .....3

    Ancillary Supply Administration Kit: Standard Syringe (Centrally Distributed).....4

    Centrally Distributed Ancillary Supply Kit Volume .....5

Direct-Ship Vaccine and Ancillary Kit.....6

    Pfizer Vaccine (Direct Ship) .....6

    Combined Ancillary Supply Kit for Administration and Mixing (Direct Ship) .....7

    Dry Ice Kit (Direct Ship).....8

    Contents and Quantities for Mega Kits to Support Pfizer Vaccine .....9

Addendum: COVID-19 Vaccine Ancillary Supply Kit Guidance ..... 11

    Possible Ancillary Kits—COVID-19 Vaccines, *Excluding Pfizer*..... 14

    Needle and Syringe Sizes for Mixing Kits..... 14

    Possible Ancillary Kits—COVID-19 Vaccine, *Pfizer Only* ..... 15

    General Needle and Syringe Resources..... 16

    Specific Ancillary Kit Products and Product Information ..... 17

## Product Information Guide Summary

### »Purpose:

This guide provides specifications for COVID-19 vaccine and associated products. **This is NOT a catalog from which you can order products.** It provides key product information including product package dimensions and weight, minimum order quantities, product presentation, distribution method, storage information, and additional information. This guide will be updated as more products become available.

### »How to use this guide

*Jurisdictional immunization program staff—use this guide to help you:*

- Learn about available COVID-19 vaccines and associated products.
- Provide answers to questions from COVID-19 vaccination providers.

*COVID-19 vaccination providers—use this guide to help you:*

- Learn about available COVID-19 vaccines and associated products.
- Prepare storage space to ensure ordered quantities can be stored at the proper temperature.
- Prepare staff responsible for receipt and storage of vaccine deliveries.

Refer to CDC’s [Vaccine Storage and Handling Toolkit](#) for recommendations and best practices on related topics such as vaccine inventory management.

### »Addendum:

The [COVID-19 Vaccine Ancillary Supply Kit Guidance](#) in the addendum (page 8) provides more specific information on the ancillary supply kits.



At this time, the products in this guide are only available for use in adults (i.e., persons aged 18 years or older). You will receive additional information, and this guide will be updated to reflect any changes COVID-19 age indications.



Additional COVID-19 vaccine may be authorized or licensed in the future. This guide will be updated as more products become available.

VTrcks: <https://www.cdc.gov/vaccines/programs/vtrcks/index.html>

Vaccine Storage and Handling Toolkit: <https://www.cdc.gov/vaccines/hcp/admin/storage/toolkit/index.html>


## Centrally Distributed Vaccines and Ancillary Kits\*

### Vaccines (Centrally Distributed)

<b>Manufacturer</b>	<a href="#">Moderna</a>
<b>NDC # in VTrckS</b>	80777-0273-99
<b>Carton dimensions</b>	53mm x 53mm x 137mm (approximately 2in x 2in x 5 3/8in)
<b>Minimum order size and increment</b>	100 doses
<b>Presentation</b>	10-dose multidose vial/10 MDV per carton
<b>Items automatically added to your VTrckS order when you select this product</b>	Ancillary supply administration kit: standard syringe
<b>Vaccination schedule</b>	2-dose series separated by at least 28 days
<b>Age indications</b>	TBD
<b>On-site vaccine storage</b>	<ul style="list-style-type: none"> <li>• -25°C to -15°C (-13°F to 5°F) in vaccine storage unit</li> <li>• 2° to 8°C (36° to 46°F) in vaccine storage unit for up to 30 days if the vial is not entered.</li> <li>• Freezer temperature settings will require adjustment if storing this vaccine and varicella-containing vaccines in the same unit. The temperature range for this vaccine is limited compared to varicella-containing vaccines.</li> </ul>

\*Pfizer's COVID-19 vaccine is **not** centrally distributed; it is distributed directly from the manufacturer. Please refer to the next section for information on Pfizer COVID-19 vaccine.

## Ancillary Supply Administration Kit: Standard Syringe (Centrally Distributed)

<b>Kit description</b>	Standard syringe kit for vaccine administration		
<b>NDC # in VTrckS</b>	11111-0001-01 (Adult)		
<b>Kit dimensions/weight</b>	14 in x 13 in x 9 in/3.5 lbs (standard ancillary adult kit)		
<b>Minimum order size and increment</b>	Kit to support administration of 100 doses (plus overage)		
<b>Accompanies 0.5mL-dose vaccines</b>	Moderna		
<b>Order Intention</b>	Initially only adult kits will be available, when authorized for use in younger populations, pediatric and mixed (pediatric and adult) will be available for ordering.		
<b>Contents: Each kit contains a label on the outside of the box with a complete inventory list.</b>	<p><b>Adult Kit</b></p> <ul style="list-style-type: none"> <li>• 85 needles (22-25G x 1")</li> <li>• 20 needles (22-25G x 1.5")</li> <li>• 105 syringes (1mL or 3mL)</li> <li>• 210 alcohol pads</li> <li>• 100 vaccination record cards</li> <li>• 1 needle gauge and length chart</li> <li>• 2 face shields</li> <li>• 4 surgical masks</li> </ul>	<p><b>Pediatric Kit</b></p> <ul style="list-style-type: none"> <li>• 105 needles (25G x 1")</li> <li>• 105 syringes (1mL or 3mL)</li> <li>• 210 alcohol pads</li> <li>• 100 vaccination record cards</li> <li>• 1 needle gauge and length chart</li> <li>• 2 face shields</li> <li>• 4 surgical masks</li> </ul>	<p><b>Mixed Kit</b></p> <ul style="list-style-type: none"> <li>• 95 needles (25G x 1")</li> <li>• 10 needles (22-25G x 1.5")</li> <li>• 105 syringes (1mL or 3mL)</li> <li>• 210 alcohol pads</li> <li>• 100 vaccination record cards</li> <li>• 1 needle gauge and length chart</li> <li>• 2 face shields</li> <li>• 4 surgical masks</li> </ul>
<b>Additional information</b>	<ul style="list-style-type: none"> <li>• Products and brands for kit components may vary.</li> <li>• All needles for vaccine administration are safety needles.</li> <li>• Due to the limited supply of needles and syringes, specification of preferences for needles or syringes is not feasible.</li> <li>• Kit and pallet configuration <ul style="list-style-type: none"> <li>○ 1 kit = 100 vaccinations (plus overage)</li> <li>○ 1 pallet = 36 kits (supports 3,600 vaccinations, plus overage)</li> </ul> </li> </ul>		
			
	<p>For more information about ancillary kits components: <a href="#">COVID-19 Vaccine Ancillary Supply Kit Guidance</a></p>		

## Centrally Distributed Ancillary Supply Kit Volume

Vaccination provider sites will need to prepare space to store vaccine and ancillary kits. Large volume orders of vaccine can be compactly packaged, however, ancillary supply kits to support these orders will require much more storage space. For each 100 doses of vaccine, an ancillary kit measuring 14 in x 13 in x 9 in and weighing 3.5 pounds will be provided. The below picture demonstrates the number and volume of ancillary kits compared to vaccine in a large volume order of 4,000 doses.

### Centrally Distributed Vaccines and Ancillary Kits



**1 Extra-large vaccine cooler  
equals  
approximately 4,000 doses  
ships with  
40 ancillary kits**



## Direct-Ship Vaccine and Ancillary Kit

### Pfizer Vaccine (Direct Ship)

<b>Manufacturer</b>	<a href="#">Pfizer</a>
<b>NDC # in VTrckS</b>	59267-1000-02
<b>Tray dimensions</b>	229mm x 229mm x 40mm (approximately 9in x 9in x 1.5in)
<b>Minimum order size and increment</b>	975 doses
<b>Presentation</b>	5-dose multidose vial/195 MDV per tray
<b>Items automatically added to your VTrckS order when you select this vaccine</b>	Combined ancillary supply kit for administration and mixing (includes 0.9% preservative-free normal saline diluent)
<b>Age indications</b>	16 years of age and older
<b>Vaccination schedule</b>	2-dose series separated by at least 21 days
<b>On-site vaccine storage</b>	<ul style="list-style-type: none"> <li>• -80°C to -60°C (-112°F and -76°F) in ultracold storage unit</li> <li>• -80°C to -60°C (-112°F and -76°F) in the original thermal shipping container with dry ice recharges. Please allow for forthcoming information on the maximum time for storage in these conditions.</li> <li>• 2° C to 8°C (36° F to 46°F) in a storage unit for up to 5 days (120 hours)</li> </ul>
<b>Additional information</b>	<ul style="list-style-type: none"> <li>• Vaccine will be shipped in a container that includes dry ice. Thermal shipping container dimensions are 400mm x 400mm x 560 mm (approximately 15 3/4in x 15 3/4in x 22in).</li> <li>• A thermal shipping container holds up to 5 cartons/trays.</li> <li>• If using the thermal shipping container to store vaccine, add dry ice pellets (9 mm to 16 mm) within 24 hours of delivery and every 5 days or as needed to maintain temperatures.</li> <li>• Unless a provider opts out, dry ice will be delivered within 24 hours of vaccine delivery to refill the thermal shipping container for the first re-ice only. Additional dry ice will not be provided. Locate a dry ice source if planning to use the shipping container to store vaccine for more than 5 days.</li> <li>• Do not use or store dry ice or liquid nitrogen (LN<sub>2</sub>) in confined areas, walk-in refrigerators, environmental chambers or rooms without ventilation. A leak in such an area could cause an oxygen-deficient atmosphere.</li> <li>• A full shipping container with vaccine and dry ice weighs approximately 80 pounds.</li> </ul>



## Combined Ancillary Supply Kit for Administration and Mixing (Direct Ship)

<b>Kit description</b>	Combined kit with small syringes for vaccine administration, mixing supplies, and diluent		
<b>NDC # in VTrckS</b>	11111-0006-02 (Adult)		
<b>Kit dimensions/weight</b>	24 in x 20 in x 24 in/40 lbs		
<b>Minimum order size and increment</b>	Kit to support administration of 975 doses (plus overage), including 0.9% preservative-free normal saline diluent		
<b>Accompanies 0.3mL-dose vaccines</b>	Pfizer		
<b>Order Intention</b>	Initially only adult kits will be available, when authorized for use in younger populations, pediatric and mixed (pediatric and adult) will be available for ordering.		
<b>Contents: Each kit contains a label on the outside of the box with a complete inventory list.</b>	<p><u>Adult Kit*</u></p> <ul style="list-style-type: none"> <li>829 needles (22-25G X 1")</li> <li>200 needles (22-25G X 1.5")</li> <li>205 mixing needles (21-25G X1.5")</li> <li>1024 syringes (1mL)</li> <li>205 syringes (3mL or 5mL)</li> <li>2458 alcohol pads</li> <li>1000 vaccination record cards</li> <li>10 needle gauge and length charts</li> <li>20 face shields</li> <li>40 surgical masks</li> <li>200 Diluent vials</li> </ul>	<p><u>Pediatric Kit*</u></p> <ul style="list-style-type: none"> <li>1024 needles (25G X 1")</li> <li>205 mixing needles (21-25G X1.5")</li> <li>1024 syringes (1mL)</li> <li>205 syringes (3mL or 5mL)</li> <li>2458 alcohol pads</li> <li>1000 vaccination record cards</li> <li>10 needle gauge and length charts</li> <li>20 face shields</li> <li>40 surgical masks</li> <li>200 Diluent vials</li> </ul>	<p><u>Mixed Kit*</u></p> <ul style="list-style-type: none"> <li>926 needles (25G X 1")</li> <li>100 needles (22-25G X1.5")</li> <li>205 mixing needles (21-25G X1.5")</li> <li>1024 syringes (1mL)</li> <li>205 syringes (3mL or 5mL)</li> <li>2458 alcohol pads</li> <li>1000 vaccination record cards</li> <li>10 needle gauge and length charts</li> <li>20 face shields</li> <li>40 surgical masks</li> <li>200 Diluent vials</li> </ul>
<b>Additional information</b>	<ul style="list-style-type: none"> <li>Products and brands for kit components may vary.</li> <li>All needles for vaccine administration are safety needles.</li> <li>Mixing needles are conventional.</li> <li>Due to the limited supply of needles and syringes, specification of preferences for needles or syringes is not feasible.</li> <li>Kit and pallet configuration               <ul style="list-style-type: none"> <li>1 kit = 975 vaccinations (plus overage)</li> <li>1 pallet = 8 kits (supports 7,800 vaccinations, plus overage)</li> </ul> </li> </ul>		
	<div style="display: flex; justify-content: space-around;">   </div>		
	<p>For more information about ancillary kit components: <a href="#">COVID-19 Vaccine Ancillary Supply Kit Guidance</a></p>		

\*The total number of needles in each adult, pediatric, and mixed combined ancillary supply kit for administration and mixing will vary (adult–1,029; pediatric–1,024; mixed–1,026). The difference in total number is related to increasing the kit build efficiency.



## Dry Ice Kit (Direct Ship)

<b>Kit description</b>	Dry ice kit with starter materials
<b>NDC # in VTrckS</b>	11111-0007-03
<b>Kit dimensions/weight</b>	TBD
<b>Minimum order size and increment</b>	Kit to support initial dry ice recharge of one thermal shipping container
<b>Accompanies ultracold vaccine</b>	Pfizer
<b>Contents</b>	<ul style="list-style-type: none"><li>• Dry ice (~35 pounds)</li><li>• Gloves for working with dry ice (1 pair)</li><li>• Face shield (1)</li><li>• Ice scoop (1)</li><li>• OSHA dry ice safety card (1)</li></ul>
<b>Additional information</b>	<ul style="list-style-type: none"><li>• Products and brands for kit components may vary.</li></ul> <p>-</p>

## Contents and Quantities for Mega Kits to Support Pfizer Vaccine

December 11, 2020

This is a clarification about the needle and syringes included in the COVID-19 ancillary supply kits for administering Pfizer’s ultra-frozen vaccine. Based on questions coming into our call centers supporting the distribution of the ancillary supply kits for Operation Warp Speed, there is some confusion about the quantities of needles and syringes included in the kits. Each kit is clearly labeled with a complete inventory list. (See example below.)



Procuring the sheer number of needles and syringes needed to support an operation of this size involves complex logistics, which is why vaccine administrators may see variance in how the supplies are packaged. The listed products may be received as a combo, defined as a single packaged needle and syringe, or as individual components. The quantities listed on the label represent the minimum number of total combinations in each kit.

- For example, if the kit above was shipped with (829) 1ml-25Gx1” syringe/needle combos and (200) 22Gx1.5” needles for administration, the kit would only contain (200) 1ml syringes to accompany the needles. This would yield a total of (1,029) needles with accompanying syringes (some packaged together, some individually).
- If a kit contained only individual needles and syringes, there may be 5 more needles than syringes. This is intentional and is designed to improve your experience and maximize the inventory across all kits built.
  - In this scenario, the adjustment was to the 22-25Gx1.5” needles. The kit was originally configured to have 195. This was increased by 5 (to 200) to improve the administrators’ experience by using two whole boxes.

Please note there will be multiple configurations used as Operation Warp Speed continues to procure products to support the assembly of these and other ancillary supply kits. Regardless of configuration, the kit contents will meet the minimum quantities listed on the kit content labels and will support the same total number of doses.

We hope this communication helps your teams as they review the content in the kits to ensure their readiness for administration when you receive the vaccine doses.

Please direct any issues regarding the ancillary kits to support the Pfizer COVID-19 vaccine to McKesson Customer Service at:

McKesson Customer Service  
Phone #: 833-272-6634  
Email: [SNSSupport@McKesson.com](mailto:SNSSupport@McKesson.com)

## **Addendum: COVID-19 Vaccine Ancillary Supply Kit Guidance**

# COVID-19 Vaccine Ancillary Supply Kit Guidance

## »Purpose:

This guidance provides an overview of the COVID-19 vaccine ancillary supply kits the U.S. Department of Health and Human Services (HHS) is providing to enrolled COVID-19 vaccination providers as part of the federal COVID-19 vaccination program. It includes a general description of seven different COVID-19 vaccine ancillary kit configurations, as well as a list of select corresponding products, product descriptions, and product quantities. For reference and to access additional product-specific information and training resources, this document also includes website hyperlinks and contact information for select product manufacturers.

## »Background:

HHS is providing ancillary supply kits for the administration of COVID-19 vaccine. The Strategic National Stockpile (SNS), managed by the HHS Office of the Assistant Secretary for Preparedness and Response (ASPR), is partnering with McKesson Corporation to produce, store and distribute these vaccine ancillary supply kits on behalf of the SNS.

## »How to use this guide:

*Jurisdictional immunization program staff—use this guide to help you:*

- Learn about products that will arrive in ancillary kit(s).

*COVID-19 vaccination providers—use this guide to help you:*

- Learn about products that will arrive in ancillary kit(s).
- Prepare storage space for the ancillary kit products you will receive.

## »Ancillary kit basics:

*How does a provider order ancillary kits?*

You do not need to order ancillary kits. When you order COVID-19 vaccine in VTrckS, ancillary supplies will automatically be ordered in amounts to match the vaccine orders.

Note: To receive COVID-19 vaccine, vaccination provider facilities and organizations must enroll in the federal COVID-19 Vaccination Program coordinated through their jurisdiction's immunization program.

*How much do ancillary kits cost?*

Ancillary kits will be provided at no cost to enrolled COVID-19 vaccination providers.

*What products are in ancillary kits?*

Refer to page 11 for possible ancillary kits.

Each ancillary kit, **except for the kits designated for use with the Pfizer vaccine**, contains supplies to administer 100 doses of vaccine (with some overage). The **Pfizer vaccination kit** will support 975 doses of vaccine (with some overage). Refer to page 12 for possible variations on kit contents.

*How will needles and syringes be packaged?*

Needles and syringes for vaccine administration may be packaged as integrated units (i.e. NO assembly required) or as separate items in a kit (i.e. assembly required).

*Can I order specific brand(s) of needles and/or syringes?*

Due to a limited supply of needles and syringes, specification of preferences for needles or syringes is not feasible. Products included in the kits may vary over time. In order to meet the demand for supplies, the federal government has purchased single-use, sterile needles and syringes from multiple manufacturers (to include foreign sources) to ensure adequate supplies. These products are approved by the Food and Drug Administration, safety-engineered and compliant with standards established by the Occupational Safety and Health Administration.

*What is considered “pediatric” and “adult”?*

For the purpose of immunizations **pediatric is birth through 18 years** and **adult is 19 years and older**. Refer to each product’s EUA for age indications.

» **Related guidance and resources:**

- [Preliminary Strategy for Distributing a COVID-19 Vaccine - PDF](#)
- [COVID-19 Vaccination Program Interim Playbook for Jurisdiction Operations - PDF](#)

» **Contact for questions:**

Vaccination providers should contact the manufacturer with questions related to proper product use. Concerns related to needle/syringe quality and performance or other contents of the vaccine ancillary supply kit should be directed to the state or local health department.

<a href="#">Becton Dickinson (BD) Worldwide</a>	<a href="https://www.bd.com/en-us/">https://www.bd.com/en-us/</a>	844-823-5433
<a href="#">Cardinal Health</a>	<a href="https://www.cardinalhealth.com/en.html">https://www.cardinalhealth.com/en.html</a>	1-800-964-5227
<a href="#">HTL Strefa</a>	<a href="https://htl-strefa.com/home-page/">https://htl-strefa.com/home-page/</a>	877-660-1900
Duopross Meditech Corp	<a href="https://www.duopross.com/">https://www.duopross.com/</a>	1-800-844-1350
<a href="#">Retractable Technologies Inc. (RTI)</a>	<a href="https://retractable.com/">https://retractable.com/</a>	888-703-1010
Marathon/Smiths Medical	N/A	941-704-7864

## Possible Ancillary Kits—COVID-19 Vaccines, *Excluding Pfizer*

The following tables list the contents of each possible ancillary kit by product and quantity. All kits are configured for 100 doses with 5% surplus.

### Needle and Syringe Sizes for *Adult* Ancillary Kits

Product	Quantity
Needle (22–25G x 1")	85
Needle (22–25G x 1.5")	20
Syringe (1mL or 3mL)	105
Alcohol Pad (sterile, individually sealed)	210
Vaccination Record Card	100
<a href="#">Needle Gauge and Length Chart</a>	1
Face Shield	2
Surgical Mask	4



### Needle and Syringe Sizes for *Pediatric* Ancillary Kits

Product	Quantity
Needle (25G x 1")	105
Syringe (1mL or 3mL)	105
Alcohol Pad (sterile, individually sealed)	210
Vaccination Record Card	100
<a href="#">Needle Gauge and Length Chart</a>	1
Face Shield	2
Surgical Mask	4



### Needle and Syringe Sizes for *Mixed* Ancillary Kits

Product	Quantity
Adult/Pediatric Needle (22–25G x 1")	85
Adult Needle (22–25G x 1.5")	20
Syringe (1mL or 3mL)	105
Alcohol Pad (sterile, individually sealed)	210
Vaccination Record Card	100
<a href="#">Needle Gauge and Length Chart</a>	1
Face Shield	2
Surgical Mask	4



## Needle and Syringe Sizes for Mixing Kits

### Needle and Syringe Sizes for *Mixing* Kits

Product	Quantity
Needle (19-21G x 1.5")	11
Syringe (5mL or 6mL)	11
Alcohol Pad (sterile, individually sealed)	22




Needle Gauge and Length Chart: [www.cdc.gov/vaccines/hcp/admin/downloads/vaccine-administration-needle-length.pdf](http://www.cdc.gov/vaccines/hcp/admin/downloads/vaccine-administration-needle-length.pdf)



## Possible Ancillary Kits—COVID-19 Vaccine, Pfizer Only

### Needle and Syringe Sizes for Pfizer Mega *Adult* Ancillary Kits

Product	Quantity
Needle (22–25G x 1")	829
Needle (22–25G x 1.5")	200
Needle, Mixing (21-25G x 1.5")	205
Syringe (1mL)	1,024
Syringe, Mixing (3mL or 5mL)	205
Alcohol Pad (sterile, individually sealed)	2,458
Vaccination Record Card	1,000
<a href="#">Needle Gauge and Length Chart</a>	10
Face Shield	20
Surgical Mask	40
Diluent	200

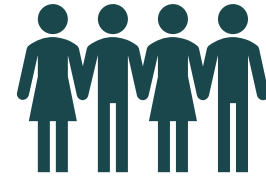

### Needle and Syringe Sizes for Pfizer Mega *Pediatric* Ancillary Kits

Product	Quantity
Needle (25G x 1")	1,024
Needle, Mixing (21-25G x 1.5")	205
Syringe (1mL)	1,204
Syringe, Mixing (3mL or 5mL)	205
Alcohol Pad (sterile, individually sealed)	2,458
Vaccination Record Card	1,000
<a href="#">Needle Gauge and Length Chart</a>	10
Face Shield	20
Surgical Mask	40
Diluent	200



### Needle and Syringe Sizes for Pfizer Mega *Mixed* Ancillary Kits

Product	Quantity
Needle (25G x 1")	926
Needle (22–25G x 1.5")	100
Needle, Mixing (21–25G x 1.5")	205
Syringe (1mL)	1,024
Syringe, Mixing (3ml or 5ml)	205
Alcohol Pads (sterile, individually sealed)	2,458
Vaccination Record Card	1,000
<a href="#">Needle Gauge and Length Chart</a>	10
Face Shield	20
Surgical Mask	40
Diluent	200



Needle Gauge and Length Chart: [www.cdc.gov/vaccines/hcp/admin/downloads/vaccine-administration-needle-length.pdf](http://www.cdc.gov/vaccines/hcp/admin/downloads/vaccine-administration-needle-length.pdf)

## General Needle and Syringe Resources



### Web Resources

- [Vaccine Administration Resource Library](#): Includes resources for preparing, administering, and documenting vaccines



### PDF Resources

- [Vaccine Administration: Needle Gauge and Length–PDF](#)
- [Vaccine Administration: Intramuscular \(IM\) Injection Children 7 through 18 years of age–PDF](#)
- [Vaccine Administration: Intramuscular \(IM\) Injection Adults 19 years of age and older–PDF](#)



### Video Resources

- [Assemble a Manufacturer-filled Syringe](#): This training addresses how to assemble a manufacturer-filled syringe, available for a variety of vaccines.

## Specific Ancillary Kit Products and Product Information

The following tables, organized by manufacturer or broker, list possible ancillary kit products with an example image (colors may vary) and links to more information. Resources for additional information are organized by:



Web resources













PDF resources



Video resources

### [Becton Dickinson Products](#) (Phone: 844-823-5433)

Item	Where to Find more Information
<p>Hypodermic Needle PrecisionGlide™ Conventional*</p> <ul style="list-style-type: none"> <li>Gauge: 22-25/Length: 1.5"</li> </ul> 	 <p><a href="#">PrecisionGlide™ Conventional Needle overview</a></p> <p>Visit <a href="#">Becton Dickinson catalog</a> and search product numbers: 305127, 305156, 305194</p>
<p>Safety Hypodermic Needle BD Eclipse™</p> <ul style="list-style-type: none"> <li>Gauge: 22-25/Length: 1-1.5"</li> <li>Gauge: 23/25/Length: 1-1.25"</li> </ul> 	 <p><a href="#">BD Eclipse™ Needle overview</a></p> <p>Visit <a href="#">Becton Dickinson catalog</a> and search product numbers: 22-25G, 1-1.5": 305762, 305761, 305763, 305767 23/25, 1-1.25": 305866, 305891, 305892</p>  <p><a href="#">BD Eclipse™ directions for use</a></p>  <p><a href="#">BD Eclipse™ Needle instruction video</a></p>
<p>Safety Combo, Syringe with BD Eclipse™ Needle 3mL BD Luer-Lok™ Syringe</p> <ul style="list-style-type: none"> <li>Gauge: 22-25/ Length: 1-1.5"</li> </ul> 	 <p><a href="#">BD Eclipse™ Needle overview</a></p> <p>Visit <a href="#">Becton Dickinson catalog</a> and search product numbers: 305782, 305783, 305787</p>  <p><a href="#">BD Eclipse™ directions for use</a></p>  <p><a href="#">BD Eclipse™ Needle Instruction video</a></p>

Item	Where to Find more Information
------	--------------------------------

1mL/3mL BD Luer-Lok™ Syringe



[Conventional syringe overview](#)

Visit [Becton Dickinson catalog](#) and search product numbers:  
309628, 309657

---

Conventional Combo\*

- Syringe (3mL/5mL) with attached needle Gauge: 22-25/Length: 1-1.5"



[Conventional syringe Overview](#)

Visit [Becton Dickinson catalog](#) and search product numbers:  
309361, 309571, 309572, 309574, 309581, 309582, 309589

---

\*For mixing ONLY and NOT for vaccine administration

PrecisionGlide™ Conventional Needle overview: <https://www.bd.com/en-us/offerings/capabilities/syringes-and-needles/conventional-syringes-and-needles/conventional-needles>

Becton Dickinson Product catalog: <http://catalog.bd.com/>

BD Eclipse™ Needle overview: <https://www.bd.com/en-us/offerings/capabilities/syringes-and-needles/safety-syringes-and-needles/safety-needles/bd-eclipse-needle>

BD Eclipse™ directions for use: [https://www.bd.com/documents/guides/directions-for-use/MPS\\_HY\\_Eclipse-needle-usage-guidelines\\_DF\\_EN.pdf](https://www.bd.com/documents/guides/directions-for-use/MPS_HY_Eclipse-needle-usage-guidelines_DF_EN.pdf)

BD Eclipse™ Needle instruction video: <https://www.bd.com/en-us/company/video-gallery?video=618677712001>

Conventional Syringe overview: <https://www.bd.com/en-us/offerings/capabilities/syringes-and-needles/conventional-syringes-and-needles/conventional-syringes>

**Cardinal Health (Phone: 1-800-964-5227)**

Item	Where to find more information
<p>Tuberculin Syringe Monoject™</p> <ul style="list-style-type: none"><li>• 1mL Luer Lock Tip</li></ul> <p>Conventional Syringe</p> 	<p></p> <p><a href="#">1 mL Luer Lock Syringes product page</a></p> <p>Visit <a href="#">Cardinal Health webpage</a> and search product numbers: 1180100777, 1180300777</p>
<p>Magellan™ 3mL Syringe with Hypodermic Safety Needle</p> <ul style="list-style-type: none"><li>• Gauge: 23/25/Length: 1"</li></ul> 	<p></p> <p><a href="#">Magellan™ Safety Needle and Syringe Combination product page</a></p> <p>Visit <a href="#">Cardinal Health webpage</a> and search product numbers: 8881833310, 8881833510</p>
<p>Magellan™ Hypodermic Safety Needle</p> <ul style="list-style-type: none"><li>• Gauge: 23/25 Length: 1"</li></ul> 	<p></p> <p><a href="#">Magellan™ Hypodermic Safety Needle Product page</a></p> <p>Visit <a href="#">Cardinal Health webpage</a> and search product numbers: 8881850310, 8881850510</p>



1 mL luer Lock Syringes product page: <https://www.cardinalhealth.com/en/product-solutions/medical/patient-care/sharp-safety/needles-and-syringes/1-ml-luer-lock-syringes.html>

Cardinal Health webpage: <https://www.cardinalhealth.com/en.html>

Magellan™ Safety Needle and Syringe Combination product page: <https://www.cardinalhealth.com/en/product-solutions/medical/patient-care/sharp-safety/needles-and-syringes/safety-needles-and-syringes/magellan-safety-needle-and-syringe-combination.html>





Magellan™ Hypodermic Safety Needle product page: <https://www.cardinalhealth.com/en/product-solutions/medical/patient-care/sharp-safety/needles-and-syringes/safety-needles-and-syringes/magellan-safety-needles.html>

**Duopross Meditech Corp (Phone: 1-800-844-1350)**

Item	Where to find more information
<p>Safety Combo, 1mL Syringe with Safety Needle*</p> <ul style="list-style-type: none"><li>Gauge: 23- 25/Length: 1-1.5"</li></ul> 	<p>Website and product information currently unavailable.</p>
<p>Safety Combo, 3mL Syringe with Safety Needle*</p> <ul style="list-style-type: none"><li>Gauge: 23-25/Length: 1-1.5"</li></ul> 	<p>Website and product information currently unavailable.</p>

\*Image for product is a placeholder stock image only. Product image is not currently available.






















**Goldbelt (Contact information currently unavailable)**

Item	Where to find more information
<p>Safety Needle*</p> <ul style="list-style-type: none"><li>• Gauge: 23-25/Length: 1-1.5"</li></ul> 	Website and product information currently unavailable.
<p>Conventional Syringe 1mL*</p> 	Website and product information currently unavailable.
<p>Conventional Syringe 3mL*</p> 	Website and product information currently unavailable.
<p>Conventional Syringe 5mL*</p> 	Website and product information currently unavailable.













\*Image for product is a placeholder stock image only. Product image is not available currently.



**Gold Coast (contact information not currently available)**


Item	Where to find more information
<p>Safety Combo, 1mL Syringe with Safety Needle</p> <ul style="list-style-type: none"> <li>Gauge: 23/Length: 1"</li> </ul> 	<div style="display: flex; justify-content: space-around;"> <div style="text-align: center;">   <a href="#">EasyTouch-FlipLock_Syringe_Instructions.pdf</a> </div> <div style="text-align: center;">   <a href="#">CP_Safety_Flyer_ProductSyringe_Needle_Brochur.pdf</a> </div> <div style="text-align: center;">   <a href="#">EasyTouch-Safety_Syringe_Instructions.pdf</a> </div> </div> <p>Product Number: Carepoint 35-8204</p>
<p>Safety Combo, 1mL Syringe with Safety Needle</p> <ul style="list-style-type: none"> <li>Gauge: 25/ Length: 1"</li> </ul> 	<div style="display: flex; justify-content: space-around;"> <div style="text-align: center;">   <a href="#">EasyTouch-FlipLock_Syringe_Instructions.pdf</a> </div> <div style="text-align: center;">   <a href="#">CP_Safety_Flyer_ProductSyringe_Needle_Brochur.pdf</a> </div> <div style="text-align: center;">   <a href="#">EasyTouch-Safety_Syringe_Instructions.pdf</a> </div> </div> <p>Product Number: Carepoint 35-8203</p>
<p>Safety Combo, FlipLock 3mL Syringe with Safety Needle</p> <ul style="list-style-type: none"> <li>Gauge: 23/Length: 1"</li> </ul> 	<div style="display: flex; justify-content: space-around;"> <div style="text-align: center;">   <a href="#">EasyTouch-FlipLock_Syringe_Instructions.pdf</a> </div> <div style="text-align: center;">   <a href="#">CP_Safety_Flyer_ProductSyringe_Needle_Brochur.pdf</a> </div> <div style="text-align: center;">   <a href="#">EasyTouch-Safety_Syringe_Instructions.pdf</a> </div> </div> <p>Product Number: 822331</p>
<p>Safety Combo, FlipLock 3mL Syringe with Safety Needle</p> <ul style="list-style-type: none"> <li>Gauge: 25/ Length: 1"</li> </ul> 	<div style="display: flex; justify-content: space-around;"> <div style="text-align: center;">   <a href="#">EasyTouch-FlipLock_Syringe_Instructions.pdf</a> </div> <div style="text-align: center;">   <a href="#">CP_Safety_Flyer_ProductSyringe_Needle_Brochur.pdf</a> </div> <div style="text-align: center;">   <a href="#">EasyTouch-Safety_Syringe_Instructions.pdf</a> </div> </div> <p>Product Number: 825231</p>
<p>Safety Needle</p> <ul style="list-style-type: none"> <li>Gauge: 25/ Length: 1"</li> </ul>  <ul style="list-style-type: none"> <li>Gauge: 23/ Length: 1"</li> </ul> 	<div style="display: flex; justify-content: space-around;"> <div style="text-align: center;">   <a href="#">EasyTouch-FlipLock_Syringe_Instructions.pdf</a> </div> <div style="text-align: center;">   <a href="#">CP_Safety_Flyer_ProductSyringe_Needle_Brochur.pdf</a> </div> <div style="text-align: center;">   <a href="#">EasyTouch-Safety_Syringe_Instructions.pdf</a> </div> </div> <p>Product Numbers:                  25G x 1": EasyTouch 802501                  23G x 1": EasyTouch 812301</p>

**HTL Strefa (Phone: 877-660-1900)**

Item	Where to find more information
<p>Safety Combo, 1mL Syringe with Safety Needle</p> <ul style="list-style-type: none"><li>Gauge: 25/Length: 1"</li></ul> 	<p> <a href="#">HTL Strefa home page</a> (product 6054 currently unlisted)</p> <p>  Appendix 1 Products picture and artwork_l:</p>
<p>Safety Combo, 3mL Syringe with Safety Needle</p> <ul style="list-style-type: none"><li>Gauge: 23/Length: 1"</li></ul> 	<p> <a href="#">HTL Strefa home page</a> (product 6053 currently unlisted)</p> <p>  Appendix 1 Products picture and artwork_l:</p>
<p>Safety Combo, 3mL Syringe with Safety Needle</p> <ul style="list-style-type: none"><li>Gauge: 25/Length: 1"</li></ul> 	<p> <a href="#">HTL Strefa home page</a> (product 6055 currently unlisted)</p> <p>  Appendix 1 Products picture and artwork_l:</p>

HTL Strefa home page: <https://htl-strefa.com/home-page/>



**Marathon/Smiths Medical (Phone: 941-704-7864)**




Item	Where to find more information
<p>Safety Combo, 3mL Syringe with Safety Needle*</p> <ul style="list-style-type: none"><li>Gauge: 22-25/Length: 1-1.5"</li></ul> 	<p>Website and product information currently unavailable.</p> <p>Product numbers: 4234, 4236, 423510</p>

\*Image for product is a placeholder stock image only. Product image is not available currently.

**Medline (Contact information currently unavailable)**

Item	Where to find more information
<p>Safety Combo, 1mL Syringe with Safety Needle</p> <ul style="list-style-type: none"><li>• Gauge: 23/Length: 1"</li></ul>	Website and product information currently unavailable.
	
<ul style="list-style-type: none"><li>• Gauge: 25/Length: 1"</li></ul>	
	



<p>Safety Needle</p> <ul style="list-style-type: none"><li>• Gauge: 23/Length: 1-1.5"</li></ul>	Website and product information currently unavailable.
	
<ul style="list-style-type: none"><li>• Gauge: 25/Length: 1"</li></ul>	
	

<p>Conventional Syringe</p> <ul style="list-style-type: none"><li>• 1mL</li></ul>	 <a href="#">Luer Lock Syringe product page</a>  Visit <a href="#">Medline webpage</a> and search product numbers: SYR101010, SYR103010
	
<ul style="list-style-type: none"><li>• 3mL</li></ul>	
	

Luer Lock Syringe product page: <https://www.medline.com/product/Luer-Lock-Syringes/Syringes-without-Needle/Z05-PF11377?question=SYR103010&index=P1&indexCount=1#mrkOrderingInfoTable>





Medline webpage: <https://www.medline.com>

**Quality Impact (Contact information currently unavailable)**

Item	Where to find more information
<p>Safety Combo, 1mL Syringe with Safety Needle*</p> <ul style="list-style-type: none"> <li>Gauge: 23-25/Size: 1"</li> </ul> 	<p>Website and product information currently unavailable.</p>
<p>Safety Combo, 3mL Syringe with Safety Needle*</p> <ul style="list-style-type: none"> <li>Gauge: 23-25/Size: 1"</li> </ul> 	<p>Website and product information currently unavailable.</p>

\*Image for product is a placeholder stock image only. Product image is not available currently.

**[Retractable Technologies Inc. \(Phone: 888-703-1010\)](#)**

Item	Where to find more information
<p>Safety Combo Unit, 1mL/3ml, with attached needle</p> <ul style="list-style-type: none"> <li>Gauge: 23/25/Length: 1"</li> </ul> 	<p>  Visit <a href="#">Retractable Technologies Inc. webpage</a> and search product numbers: 10161, 10311, 10391</p> <p>  <a href="#">VanishPoint® Syringes brochure</a>  <a href="#">VanishPoint® Syringes product usage information</a></p> <p>  <a href="#">VanishPoint® Syringes video</a></p>

Retractable Technologies Inc. webpage: <https://vanishpoint.com/>

VanishPoint® Syringes brochure: <https://d2ghdaxqb194v2.cloudfront.net/577/166728.pdf>

VanishPoint® Syringes product usage information: <https://d2ghdaxqb194v2.cloudfront.net/577/166715.pdf>

VanishPoint® Syringes video: <https://youtu.be/wC-uXq3uUdQ>